

FREEZE TIME

REVERSE AGING




VENUSFREEZE *Plus*



VENUS FREEZE *Plus*

The easiest-to-use system with advanced safety features for non-invasive wrinkle reduction via skin tightening. Treatments are pain-free with no downtime and for all skin types. It showcases two new features for enhanced patient safety, improved ease of use for the operator, and more consistent and predictable results.

NON-INVASIVE ANTI-AGING TREATMENTS

Tighten lax skin and reduce wrinkles and rhytides through neocollagenesis and fibroblast stimulation.

ADVANCED SAFETY FEATURES

Integrated Automatic Temperature Control (ATC) delivers optimal patient safety and comfort throughout the treatment.

EASIEST TO USE

Preset treatment parameters, ergonomic handpieces, and Real-Time Thermal Feedback deliver unparalleled ease of use for the operator.

NO PAIN, NO DOWNTIME AND SAFE FOR ALL SKIN TYPES

Radio Frequency (RF) technology is proven to be safe for all skin types. Together with being painless and having no downtime, you will see higher patient compliance with Venus Freeze Plus™ treatments.

HIGHER RETURN ON INVESTMENT

Venus Concept's unique business model ensures lower-risk and higher-return ownership.

APPLICATORS & (MP)² TECHNOLOGY

DIAMONDPOLAR™ AND OCTIPOLAR™

Venus Freeze Plus™ features two applicators: DiamondPolar™ for smaller treatment areas and OctiPolar™ for larger areas. Both of these lightweight, ergonomically designed applicators are powered by our proprietary (MP)² technology, combining Multi-Polar Radio Frequency and Pulsed Electro Magnetic Fields.



MULTI-POLAR RADIO FREQUENCY

Multi-Polar Radio Frequency (RF) uses a complex algorithm to deliver more homogeneous heating to multiple tissue depths, allowing for faster buildup of heat and easier maintenance of the therapeutic temperature. This results in greater efficacy and enhanced patient safety. With no heat spikes, there is also no need for topical cooling agents and patients report an enjoyable, pain-free treatment experience.

PULSED ELECTRO MAGNETIC FIELDS

The effect of Multi-Polar RF is enhanced by Pulsed Electro Magnetic Fields (PEMF), a non-thermal technology emitted through the applicator's electrodes. PEMF has been proven effective in triggering regenerative processes in the skin and has been used in conventional medicine for decades.

THE SYNERGISTIC EFFECT OF MULTI-POLAR RF AND PEMF

(MP)² synergistically increases collagen production through two independent mechanisms, Multi-Polar RF (thermal) and PEMF (non-thermal). The RF directly stimulates fibroblasts, while PEMF induces fibroblast proliferation through the release of the growth factor FGF-2. Together, they work to increase collagen synthesis. (MP)² is proven and effective in remodeling of collagen for skin tightening and the creation of new capillaries, which renews blood supply.





WHY VENUS FREEZE *Plus* IS THE

BEST OPTION FOR YOU

Venus Freeze Plus™ is the easiest to use, which makes it highly delegable. It includes two integrated safety features for enhanced patient comfort, improved ease of use for the operator, and more consistent and predictable outcomes.



REAL-TIME THERMAL FEEDBACK

Venus Freeze Plus™ applicators are equipped with integrated sensors that provide real-time thermal feedback, allowing for easy and immediate monitoring of skin temperature profile. The system's display shows the temperature profile throughout the treatment, so the operator can ensure optimal safety and effectiveness.


AUTOMATIC TEMPERATURE CONTROL

Venus Freeze Plus™ features Automatic Temperature Control (ATC), which allows the operator to easily select a target temperature suitable for the specific patient and treatment area. This unique feature makes it easier to maintain a steady therapeutic temperature throughout the treatment, enhancing the consistency and predictability of results. ATC also contributes to the safety of the treatment.



EXPERIENCE

THE RESULTS

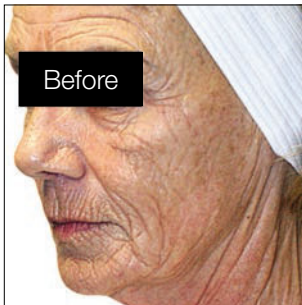


Courtesy of Dima Ali, MD
Results after 7 treatments





Courtesy of Eric Adler, MD
Results after 10 treatments



Courtesy of Scott Callahan, MD
Results after 6 treatments

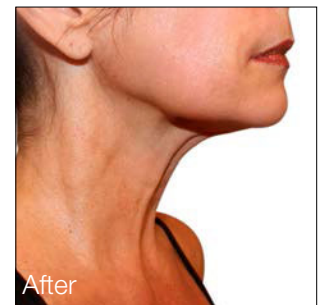

Courtesy of Scott Callahan, MD
Results after 6 treatments

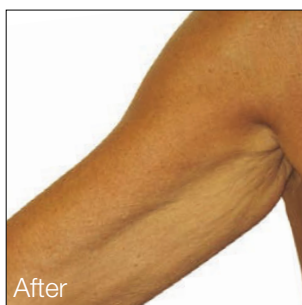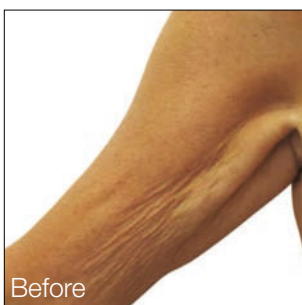

Courtesy of Anti-Aging Ltd.
Results after 6 treatments



Courtesy of Lowcountry Plastic Surgery
Results after 3 treatments

DEVICE SPECIFICATIONS

TECHNICAL SPECIFICATIONS

| | |
|--|--|
| Input Voltage | 100-240 V, 50-60 Hz, Max. 5A |
| Output RF Frequency | 1 MHz |
| Max. RF Output Power | Up to 150 watts |
| (MP) ² Applicators | DiamondPolar™ Applicator OctiPolar™ Applicator |
| Number of (MP) ² Synthesizers | Up to 8 |
| Magnetic Pulse Frequency | 15 Hz |
| Magnetic Field | 15 Gauss |
| Dimensions | 65 cm x 55 cm x 110 cm / 25.6 in x 21.7 in x 43.3 in (D x W x H) |
| Weight | 25 kg / 55 lb |



Corporate Office: 255 Consumers Road, Suite 110 | Toronto, Ontario, M2J 1R4

USA Office: 4556 N Hiatus Road | Sunrise, Florida, 33351

Tel: 888-907-0115 Visit us at venusconcept.com



Venus Freeze Plus™ is licensed by Health Canada for non-invasive skin tightening. © All rights reserved. Venus Freeze Plus™ is a trademark of Venus Concept™ and may be registered in certain jurisdictions. (MP)² is a proprietary technology of Venus Concept™.